

TAS[®] (Toe Alignment Splint)

Toe Alignment Splint



High wearing comfort,
even under socks

Elastic
heel strap

Elastic
metatarsal strap

- Universal size 
- Right = Left 
- Low weight 
- Latex free 




The TAS® toe alignment splint is a corrective brace. It is used conservatively and postoperatively in place of taping. The hallux valgus toe orthosis is an effective and inexpensive way of immobilising the operated/corrected toe in the desired position and holding it in place during the healing process. The elastic heel strap on the back of the foot prevents the midfoot splint from slipping. The individual loops hold the toes in the desired position. The soft material ensures a high level of comfort. It is also ideal for use after minimally invasive surgery.

Indications

For conservative and postoperative correction / realignment of toes in:

- > Hallux valgus (1)
- > Tailor's bunio (2)
- > Hammer toes (3)

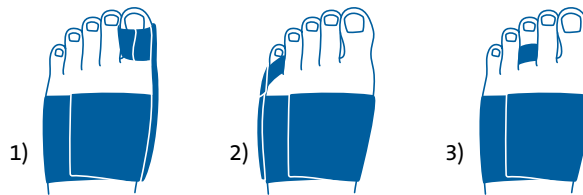

Features

- > Can be used conservatively and postoperatively
- > Redressing bandage instead of tape bandage
- > Freely positionable fixation + alignment of operated / corrected toes
- > Also for use after minimally invasive procedures
- > Soft, skin-friendly material → Pleasant wearing comfort
- > Stabilising element: Thermoplastically shapeable
- > Elastic, 34 cm long metatarsal strap → Care of wide feet possible
- > Elastic heel strap on the rearfoot
- > Can be worn with post-op and healing shoes
- > Universal fit for right and left



Adduction



Abduction



Dorsal extension



Plantar flexion

SIZE

ART.-NO.



Colour: White

Universal TAS-22



Product Video



DARCO (Europe) GmbH

Gewerbegebiet 18 | 82399 Raisting | Germany | Tel. +49 8807 9228-0 | Fax -22 | info@darco.de | www.darco.de/en

